


Super Arrow-Flex[®] Sheath Introducer

FLEXIBILITY YOU CAN SEE. STRENGTH YOU CAN FEEL.

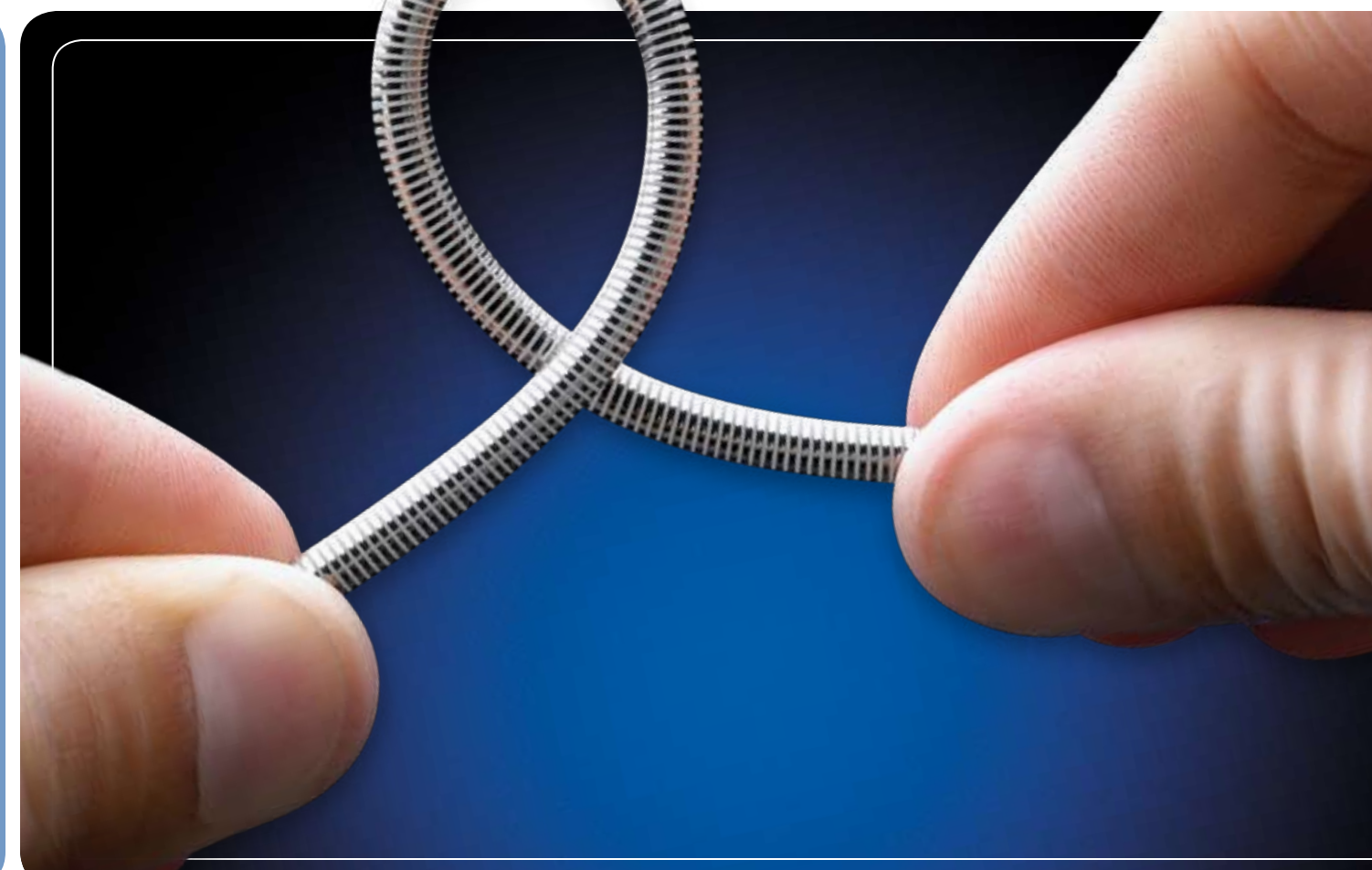

To locate an Arrow representative or for ordering information, contact Customer Service at:

United States International
Tel: 800.523.8446 Tel: +1.610.655.8522
Fax: 800.343.2935 Fax: +1.610.655.8566

www.arrowintl.com

Super Arrow-Flex Sheath Introducer Sets With Wire Guides

Product No.	Sheath Size	Sheath Length	Color Coded Hub	Sets Per Case
CP-07511-P	5 Fr.	7.5cm	Gray	10
CP-07511	5 Fr.	11cm	Gray	10
CP-07611-P	6 Fr.	7.5cm	Green	10
CP-07611	6 Fr.	11cm	Green	10
CP-07711	7 Fr.	11cm	Orange	10
CP-07811	8 Fr.	11cm	Blue	10
CP-07911	9 Fr.	11cm	Black	10
CP-07011	10 Fr.	11cm	White	10

Super Arrow-Flex Sheath Introducer Sets Without Wire Guides

Product No.	Sheath Size	Sheath Length	Color Coded Hub	Sets Per Case
5 Fr. Introducers				
CL-07511	5 Fr.	11cm		10
CL-07524	5 Fr.	24cm		5
CL-07545	5 Fr.	45cm	Gray	5
CL-07565	5 Fr.	65cm		5
CL-07590	5 Fr.	90cm		2
6 Fr. Introducers				
CL-07611	6 Fr.	11cm		10
CL-07624	6 Fr.	24cm		5
CL-07635	6 Fr.	35cm		5
CL-07645	6 Fr.	45cm	Green	5
CR-07645*	6 Fr.	45cm (Renal)		5
CL-07665	6 Fr.	65cm		5
CL-07690	6 Fr.	90cm		2
7 Fr. Introducers				
CL-07711	7 Fr.	11cm		10
CL-07724	7 Fr.	24cm		5
CL-07735	7 Fr.	35cm		5
CL-07745	7 Fr.	45cm		5
CR-07745*	7 Fr.	45cm (Renal)		5
CR-07745-NT [†]	7 Fr.	45cm (Renal)	Orange	5
CL-07765	7 Fr.	65cm		5
CL-07780	7 Fr.	80cm		2
CL-07790-R*	7 Fr.	90cm (Carotid)		2
CL-07700	7 Fr.	100cm		2
8 Fr. Introducers				
CL-07811	8 Fr.	11cm		10
CL-07824	8 Fr.	24cm		5
CL-07835	8 Fr.	35cm		5
CL-07845	8 Fr.	45cm		5
CT-07860	8 Fr.	60cm (Transseptal)	Blue	2
CL-07865	8 Fr.	65cm		5
CL-07880	8 Fr.	80cm		2
CL-07890-R**	8 Fr.	90cm (Carotid)		2
CL-07800	8 Fr.	100cm		2

Each "CP" set includes:

- One: Super Arrow-Flex Sheath with Integral Side Port/Hemostasis Valve and attached 3-Way Stopcock
- One: 17 3/4" (45cm) long dual purpose "J" tip spring-wire guide in holder with tip straightener. Short-length sheath set (P) has .021" dia. and all others have .035" dia. wire guide
- One: Vessel Dilator with Snap-Lock Feature
- One: Obturator Cap

Product No.	Sheath Size	Sheath Length	Color Coded Hub	Sets Per Case
9 Fr. Introducers				
CL-07911	9 Fr.	11cm		10
CL-07924	9 Fr.	24cm		5
CL-07965	9 Fr.	65cm	Black	5
CL-07980	9 Fr.	80cm		2
CL-07900	9 Fr.	100cm		2
10 Fr. Introducers				
CL-07011	10 Fr.	11cm		10
CL-07024	10 Fr.	24cm		5
CL-07035	10 Fr.	35cm	White	5
CL-07045	10 Fr.	45cm		5
CL-07065	10 Fr.	65cm		5
CL-07080	10 Fr.	80cm		2
11 Fr. Introducers				
CL-71165	11 Fr.	65cm	Yellow	5
CL-71180	11 Fr.	80cm		2

Each "CL" set includes:

- One: Super Arrow-Flex Sheath with Integral Side Port/Hemostasis Valve and attached 3-Way Stopcock
- One: Vessel Dilator with Snap-Lock Feature
- One: Obturator Cap

*These products include a special Tuohy-Borst "Y" Adapter with Side Arm Extension.
[†]Tuohy-Borst "Y" Adapter is not included with this product.
[‡]Contact your Arrow representative for availability.

Caution: U.S. Federal law limits this device to sale by or on order of a physician.

U.S. Patent Nos. 5,180,376 and 5,484,425.



ARROW
INTERNATIONAL

Distribution Worldwide:

Arrow offices are located in Belgium, Canada, Czech Republic, France, Germany, Greece, Hungary, India, Italy, Japan, Mexico, Netherlands, Slovakia, South Africa, Spain, United Kingdom, and the United States.

Arrow International, Inc., 2400 Bernville Road, Reading, PA 19605 U.S.A. Tel: 610.378.0131

Authorized EC representative: Arrow Deutschland GmbH, Justus-von-Liebig-Strasse 2, D-85435 Erding, Germany


www.arrowintl.com



SHEATH INTRODUCERS

ARROW
INTERNATIONAL

Complete vascular access without kinking or collapsing

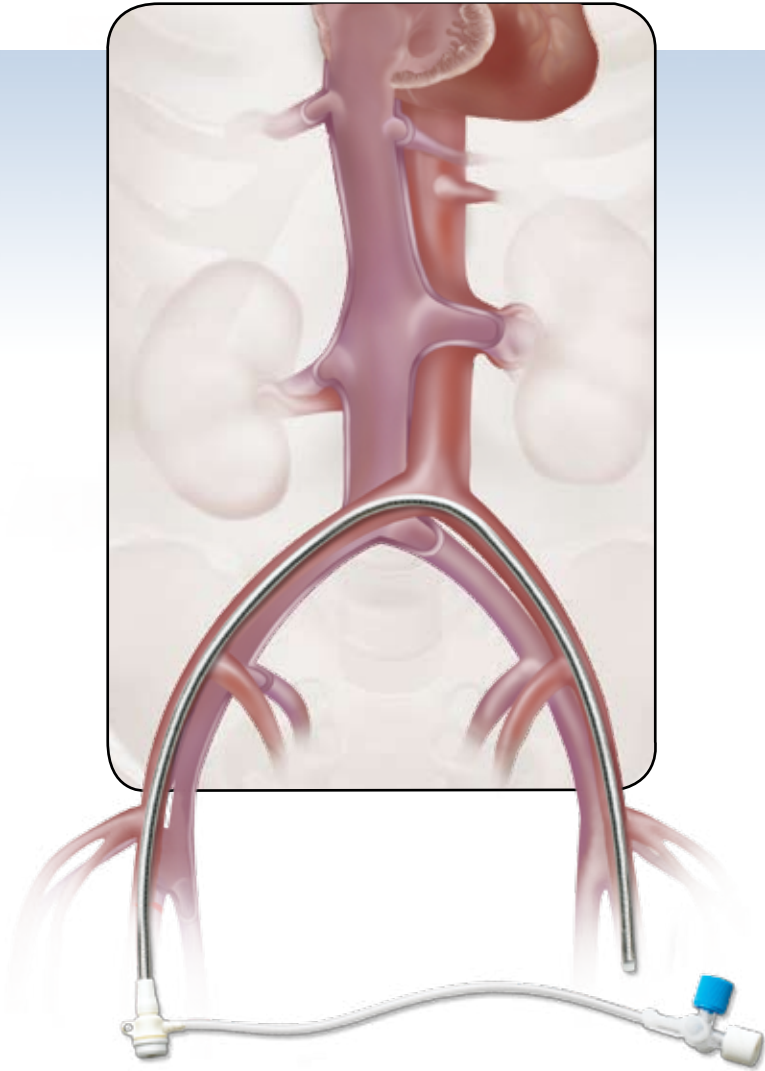

THE FLEXIBLE, KINK-RESISTANT SUPER ARROW-FLEX SHEATH INTRODUCER

Super Arrow-Flex gives you complete vascular access, thanks to its superior design and construction. From its easy insertion and handling to its highly visible radiopaque tip,* Super Arrow-Flex is designed to deliver optimal performance from start to finish. And with a wide range of sheath sizes and lengths, Super Arrow-Flex gives you flexibility in choices for multiple vascular applications.

MAXIMUM EFFECTIVENESS EVEN IN THE MOST CHALLENGING CASES

- Exclusive coil-wire design allows the sheath to flex at any point in any direction without kinking or losing support
- Excellent steerability helps negotiate tortuous anatomies
- Tapered dilator provides exceptional pushability to the contralateral leg and the Superficial Femoral Artery (SFA)



EASY PLACEMENT AND INSERTION

- Radiopaque tip marker enhances visualization during insertion and enables precise deployment of a catheter just beyond the sheath tip*
- Highly radiopaque sheath, yet spacing between coils enables visualization through the sheath
- Smooth transition between the sheath and dilator helps to prevent peel-back
- Hydrophilic surface helps facilitate smooth insertion
- Valve lubrication eases catheter placement through the sheath



RANGE OF PRODUCTS FOR INTERVENTIONAL PROCEDURES

- Specific designs suitable for renal, carotid, and transeptal procedures
- Wide range of products providing reliable access for crossover, TIPS, and antegrade brachial procedures
- Variety of product options from 5 Fr. to 11 Fr. and 7.5 cm to 100 cm
- Sheath available with and without guide wire
- Color coded for convenience

SAFETY AND PEACE OF MIND

- Ability to be used as a guiding sheath may eliminate the need for a guiding catheter (includes obturator cap)
- Color-coded hubs for easy identification of sheath sizes
- Tri-cuspid hemostasis valve facilitates catheter manipulation with minimal bleed-back
- Sheath and dilator lock together to enable simultaneous, secure insertion



*In all 35 cm and longer except 5 Fr.