


# Perfecting Performance



**Occlutech® Delivery Set (ODS)** Perfect Delivery



# Occlutech® Delivery Set (ODS)

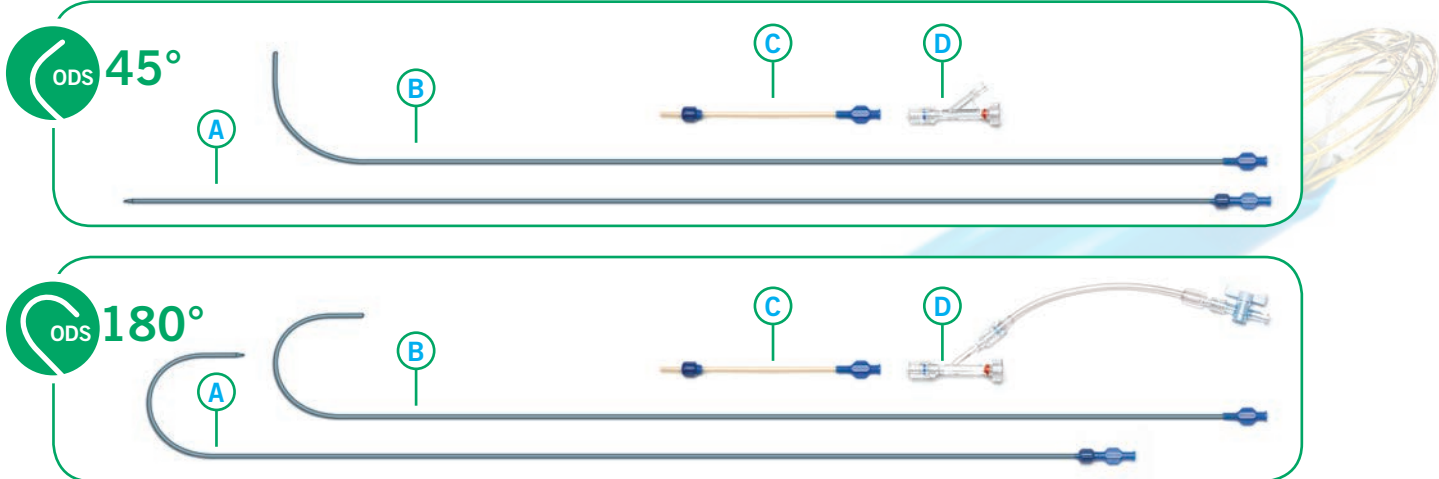

## Features & Benefits

- A. Dilator** – The dilator (A) is made from smooth conformable low friction radio-opaque material. The dilator is inserted into the sheath (B) to form a supporting inner body. The inner lumen of the dilator has 0,035" guide wire compatibility.
- B. Delivery sheath** – The sheath (B) is of Mullins type and consists of a smooth outer layer and a low friction inner layer reinforced with metal braiding. The braid improves the strength of the sheath and provides radiopacity. The sheath has a proximal female connector for connection to the dilator (A) or to the loader (C).
- C. Transparent loader** – The loader (C) is made of low friction material and transparent to allow inspection for air bubbles after pulling the occluder into the loader.
- D. Y-connector with hemostasis valve** – The Y-connector (D) is a commercially available connector with a side-port and hemostasis valve. The valve allows for leak-proof passage of the Occluder delivery wire. The side-port permits flushing before and during the implant procedure.



## Ordering information – dimensions and sizes

### Occlutech® Delivery Set 45°

Article No.	Size (F)	Length Sheath (mm)	Length Loader (mm)	Dilator Tip (mm)	Hemo-valve (mm)	Distal Angle
51DS007	7	800	120	25	70	45°
51DS008	8	800	160	25	70	45°
51DS009	9	800	160	25	70	45°
51DS010	10	800	160	25	70	45°
51DS011	11	800	200	25	70	45°
51DS012	12	800	200	25	70	45°
51DS014	14	800	200	25	70	45°

### Occlutech® Delivery Set 180°

Article No.	Size (F)	Length Sheath (mm)	Length Loader (mm)	Dilator Tip (mm)	Hemo-valve (mm)	Distal Angle
51US006	6	800	120	25	70	180°
51US007	7	800	120	25	70	180°
51US008	8	800	160	25	70	180°
51US009	9	800	160	25	70	180°

**Customer Service International**  
 Occlutech, La Cours Gata 2, S-252 31 Helsingborg, Sweden  
 Tel +46 (0)42 33 65 32, Fax +46 (0)42 311 09 70  
 order@occlutech.com, [www.occlutech.com](http://www.occlutech.com)