


Vascular Intervention // Peripheral

Passeo-35 HP

High Pressure PTA Balloon/0.035"/OTW

Indicated for arterio-venous dialysis fistulae

Treatment of femoral, iliac and renal arteries

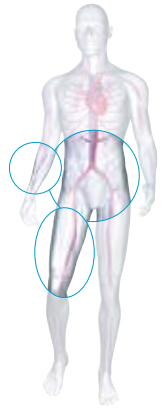

- Up to 27 atm to dilate resistant lesions in complex anatomy
- 96 % success rate opening arterio-venous dialysis fistulae¹
- Controlled compliance for excellent stability during inflation
- Flexible and conformable minimizing vessel straightening

¹ BIOTRONIK data on file (PME report, 2013)

Passeo-35 HP

High Pressure PTA Balloon/0.035"/OTW

Technical Data	Balloon catheter
Catheter type	OTW
Recommended guide wire	0.035"
Tip	Soft, short, tapered
Balloon material	Nylon/Pebax, controlled compliance
Balloon folding	3-fold (ø 3.0 - 9.0 mm); 5-fold (ø 10.0 - 12.0 mm)
Balloon markers	2 radiopaque markers
Sizes	ø 3.0 - 12.0 mm; L: 20 - 100 mm
Shaft	5.9F, coaxial
Usable length	40 and 75 cm

Compliance Chart		Balloon diameter x length (mm)								
		ø 3.0 x 40	ø 4.0 x 20-40	ø 5.0 x 20-60	ø 6.0 x 20-100	ø 7.0 x 20-100	ø 8.0 x 20-80	ø 9.0 x 40	ø 10.0 x 40	ø 12.0 x 40
Nominal Pressure (NP)	atm*	14	14	14	14	14	14	12	12	12
	ø (mm)	3.11	4.01	5.01	6.05	6.93	7.98	8.96	10.02	11.86
Rated Burst Pressure (RBP)	atm*	27	27	27	25	23	22	20	20	18
	ø (mm)	3.42	4.41	5.46	6.56	7.45	8.50	9.66	10.78	12.41

* 1 atm = 1.013 bar

Ordering Information	Balloon ø (mm)	Catheter length 40 cm					Catheter length 75 cm				
		Balloon length (mm)					Balloon length (mm)				
		20	40	60	80	100	20	40	60	80	100
	3.0	-	399056	-	-	-	-	399077	-	-	-
	4.0	399057	399058	-	-	-	399078	399079	-	-	-
6F	5.0	399059	399060	399061	-	-	399080	399081	399082	-	-
	6.0	399062	399063	399064	-	399065	399083	399084	399085	-	399086
	7.0	399066	399067	399068	-	399069	399087	399088	399089	-	399090
	8.0	399070	399071	399072	399073	-	399091	399092	399093	399094	-
7F	9.0	-	399074	-	-	-	-	399095	-	-	-
	10.0	-	399075	-	-	-	-	399096	-	-	-
8F	12.0	-	399076	-	-	-	-	399097	-	-	-

Passeo-35 HP is part of the BIOTRONIK **6F** Solutions portfolio, including:

■ Introducer Sheath: **Fortress** ■ Balloons: **Passeo-35, AngioSculpt** ■ Stents: **Dynamic, Astron, Pulsar-35**

Manufacturer:
Creagh Medical
IDA Business Park
Ballinasloe
Co. Galway
Ireland

Distributor:
BIOTRONIK AG
Ackerstrasse 6
8180 Bülach · Switzerland
Tel +41 (0)44 8645111
Fax +41 (0)44 8645005
info.vi@biotronik.com
www.biotronik.com

For ordering please contact your
local sales representative