


# Passeo-35

PTA Balloon/0.035"/OTW

Indicated for lower limb, renal<sup>1</sup> and iliac<sup>1</sup> arteries

Treatment of arterio-venous dialysis fistulae

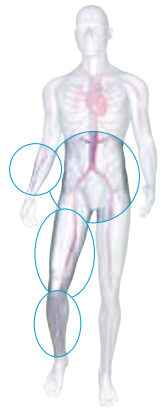

- Low profile, wide range:  $\varnothing$  3 - 10 mm and up to 200 mm balloon length
- Puncture resistant balloon material for most calcified lesions
- 5F shaft for uncompromised flush capability and angiographic quality
- Patchwork coating for enhanced crossability

<sup>1</sup> Australia: not TGA approved for use within the renal and common iliac arteries.

# Passeo-35 – PTA Balloon/0.035"/OTW

## Technical Data

Balloon catheter	
Catheter type	OTW
Recommended guide wire	0.035"
Tip	Soft, short, tapered, colored
Balloon material	SCP (Semi Crystalline Polymer), controlled compliance
Balloon folding	5-fold
Balloon coating	Hydrophobic patchwork coating
Balloon markers	2 swaged markers
Sizes	ø 3.0 - 10.0 mm; L: 20 - 200 mm
Shaft	5F, hydrophobic coating, dual-lumen
Usable length	80, 90 and 130 cm
Guide wire lumen	Hydrophobic coating

## Compliance Chart

		Balloon diameter x length (mm)													
		ø 3.0 x 20-40	ø 3.0 x 60-200	ø 4.0 x 20-40	ø 4.0 x 60-200	ø 5.0 x 20-100	ø 5.0 x 120-200	ø 6.0 x 20-100	ø 6.0 x 120-200	ø 7.0 x 20-100	ø 7.0 x 120-200	ø 8.0 x 20-100	ø 9.0 x 20-80	ø 10.0 x 20-40 <sup>3</sup>	ø 10.0 x 20-80
Nominal Pressure (NP)	atm <sup>2</sup>	7	7	7	7	7	7	7	7	7	7	7	7	7	7
	ø (mm)	3.0	3.00	4.0	4.00	5.0	5.0	6.0	6.0	7.0	7.0	8.0	9.0	10.0	10.0
Rated Burst Pressure (RBP)	atm <sup>2</sup>	20	20	20	18	16	16	16	16	14	14	14	12	12	11
	ø (mm)	3.5	3.20	4.5	4.32	5.7	5.3	6.7	6.4	7.7	7.3	8.7	9.4	10.5	10.3

<sup>3</sup> Usable length: 80 cm; <sup>2</sup> 1 atm = 1.013 bar

## Ordering Information

Balloon ø (mm)	Catheter length 80 cm Balloon length (mm)					Catheter length 90 cm Balloon length (mm)							
	20	40	60	80	100	60	80	100	120	150	170	200	
5F	3.0	359545	359547	-	-	-	383231	383235	383239	383243	389775	389776	387162
	4.0	359546	359548	-	-	-	383232	383236	383240	383244	383248	383252	383256
	5.0	357282	357288	357294	357298	357302	-	-	-	383245	383249	383253	383257
	6.0	357283	357289	357295	357299	357303	-	-	-	383246	383250	383254	383258
	7.0	357284	357290	357296	357300	357304	-	-	-	383247	383251	383255	383259
6F	8.0	357285	357291	357297	357301	357305	-	-	-	-	-	-	-
	9.0	357286	357292	-	-	-	383233	383237	-	-	-	-	-
	10.0	357287	357293	-	-	-	383234	383238	-	-	-	-	-

Balloon ø (mm)	Catheter length 130 cm Balloon length (mm)									
	20	40	60	80	100	120	150	170	200	
5F	3.0	359549	359551	383264	383268	383272	383276	389777	389778	387163
	4.0	359550	359552	383265	383269	383273	383277	383281	383285	383289
	5.0	357306	357310	357314	357318	357322	383278	383282	383286	383290
	6.0	357307	357311	357315	357319	357323	383279	383283	383287	383291
	7.0	357308	357312	357316	357320	357324	383280	383284	383288	383292
6F	8.0	357309	357313	357317	357321	357325	-	-	-	-
	9.0	383260	383262	383266	383270	-	-	-	-	-
	10.0	383261	383263	383267	383271	-	-	-	-	-

Passeo-35 is part of the BIOTRONIK **6F** Solutions portfolio, including:

■ Introducer Sheath: **Fortress** ■ Balloons: **Passeo-35 HP, AngioSculpt** ■ Stents: **Dynamic, Astron, Pulsar-35**

For ordering please contact your  
local sales representative

BIOTRONIK AG  
Ackerstrasse 6  
8180 Bülach · Switzerland  
Tel +41 (0)44 8645111  
Fax +41 (0)44 8645005  
info.vi@biotronik.com  
www.biotronik.com