



## High Torque

### Shaft braiding



#### Conventional catheter

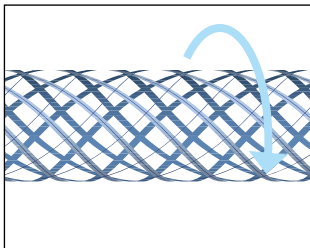

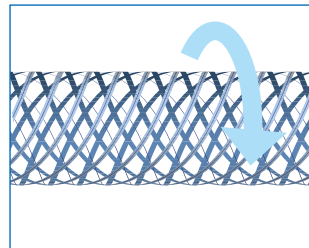




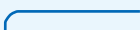
When the braiding contained within the catheter body is positioned horizontally, the torquability of the catheter is reduced.



#### HT (High Torque) catheter

By ensuring the braiding is positioned vertically, the torquability of the catheter is maximized.

### Ordering information

Product Type	Shape	Length	Side Hole	Model Number		
				4Fr	5Fr	
Judkins Right		JR3.5	100	0	-	HT5-JR35A
		JR4.0	100	0	-	HT5-JR40A
		JR5.0	100	0	-	HT5-JR50A
Judkins Left		JL3.5	100	0	-	HT5-JL35A
		JL4.0	100	0	-	HT5-JL40A
		JL5.0	100	0	-	HT5-JL50A
Radial		Radial 3.5	100	1	HT4-LRT35	HT5-LRT35
		Radial 4.0	100	2	-	HT5-LRT40

Manufacturer



<http://www.goodmankk.com/english/>

- ▶ Please consult your sales representative regarding stock availability.
- ▶ Goodtec is a trademark of Goodman Co., Ltd.
- ▶ Please refer to the Instructions For Use for further details regarding product use and specifications.