



Profit[™] 6F Plus
Guiding Catheter

Gain the advantage

 GOODMAN

Profit™ 6F Plus

Guiding Catheter

Features

0.0072 inch Inner Diameter

A large internal diameter (0.0072 inch) and PTFE liner provides a low frictional coefficient, making it ideal for large devices and performing KBT technique.

Optimal backup force

A gradual transition in rigidity from tip to shaft ensures gentle ostial engagement while maintaining strong back-up force.


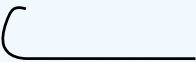



Superior control

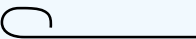

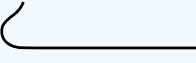
While offering strong kink resistance and remarkable operability, the flexible shaft provides for required support during the procedure.

Safe soft tip

Along with a flexible and radiopaque tip, a 'rounding' process increases the safety profile when engaging the vessel.

Ordering information

Product Type	Shape	Model Number
Judkins Right 	JR	P6P-JR35
		P6P-JR40
		P6P-JR35SH
		P6P-JR40SH
Radial Unlimited 	RU	P6P-RU
		P6P-RUM375
		P6P-RUSH
Radial Right 	RR	P6P-RR
Short Amplatz Left 	SAL	P6P-SAL1
		P6P-SAL1SH
Judkins Left 	JL	P6P-JL35
		P6P-JL40
		P6P-JL45
		P6P-JL35SH
		P6P-JL40SH
		P6P-JL40SH

Product Type	Shape	Model Number
Judkins Left Short Tip 	JL	P6P-JL35ST
		P6P-JL40ST
Special Support 	SS	P6P-SS325
		P6P-SS35
		P6P-SS375
		P6P-SS40
		P6P-SS325SH
Amplatz Left 	AL	P6P-SS35SH
		P6P-AL075
		P6P-AL1
		P6P-AL15
		P6P-AL2
		P6P-AL075SH
P6P-AL1SH		

- ▶ Please consult your sales representative regarding stock availability.
- ▶ Please refer to the Instructions For Use for further details regarding product use and specifications.