


Conical Model

Conical Model

32 x 26	150 FF 25
34 x 28	150 FF 25
36 x 30	150 FF 25
38 x 32	150 FF 25
40 x 34	150 FF 25
42 x 36	150 FF 25
44 x 38	150 FF 25
46 x 40	150 FF 25

Free Flow - FF25 (With coated protection)

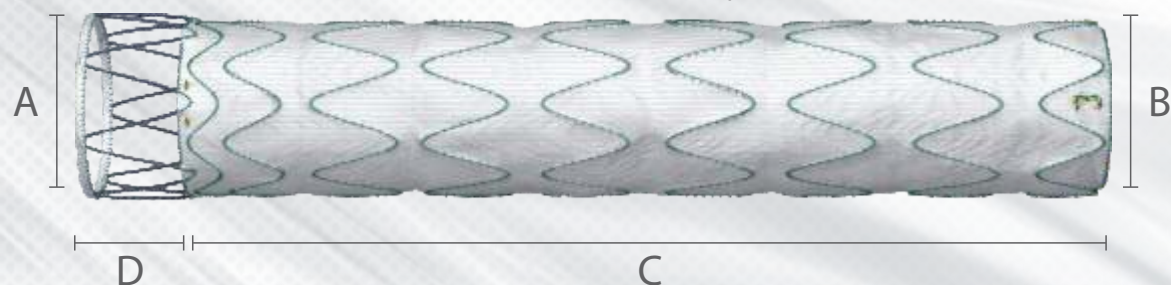

Straight Model

A	B	C	D
Proximal Diameter 28 a 46 mm	Distal Diameter 26 a 46 mm	Length 110 - 130 - 150 - 190 mm	Finishing Free Flow 25 mm
Profile 18 e 20 Fr			


***See availability of customization/other measures**



AGENCIAPROSPECTA.COM

DOMINUS

Endoprótese *Stent-Graft* Dominus

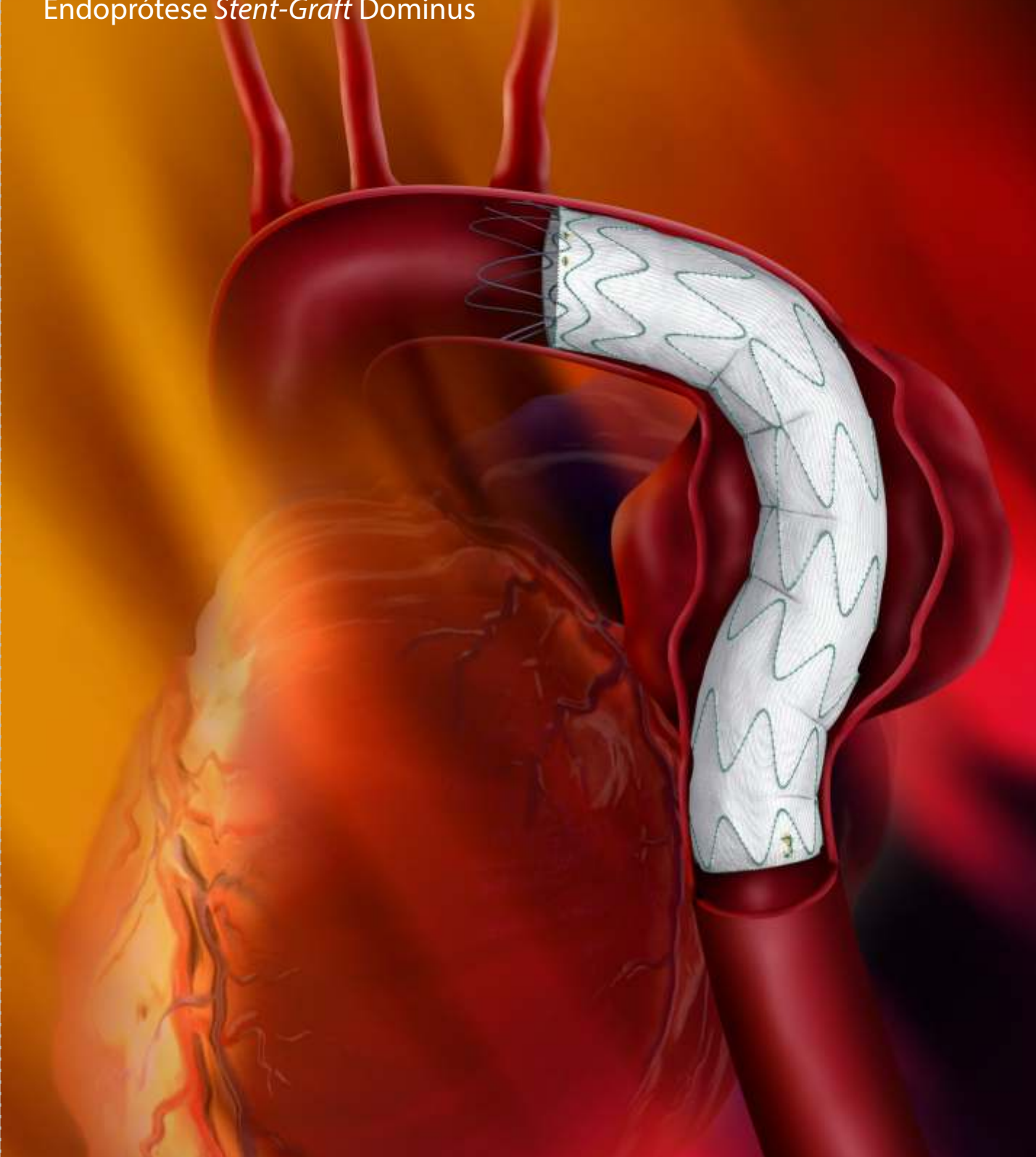

Braile Biomédica Indústria, Comércio e Representações Ltda.

Av. Pres. Juscelino K. de Oliveira, 1505 • Jardim Tarraf I
CEP 15091-450 • São José do Rio Preto • SP • Brasil
Telefone 17 2136 7000 • Fax 17 2136 7040
Other countries • Phone # (55) 17 2136 7000 • Fax 17 2136 7040
CNPJ: 52.828.936/0001-09
Responsável Técnico: Vladimir D. A. Ramirez • CRF/SP 9010
Registro ANVISA nº 10159030093
www.braile.com.br

BRAILE
BIOMÉDICA


www.braile.com.br

BRAILE
BIOMÉDICA




► **Free Flow**

Non-coated stent for expansion of attaching area

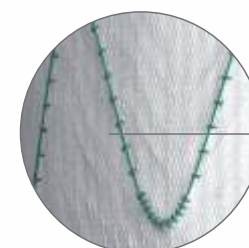

► **Suture**

Sutures made of polyester for attaching the metallic framework

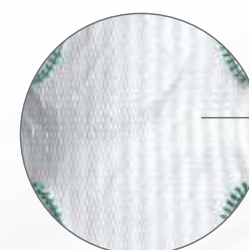


► **Coating**

Low porosity tubular polyester

► **Nitinol (nickel titanium)**

Superelastic metal alloy with thermal memory

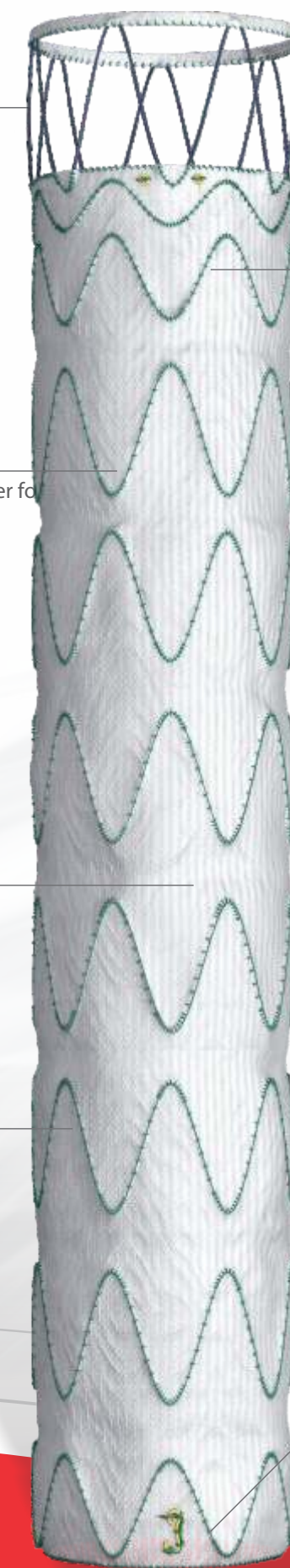


► **Radial Force**

Proximal and distal areas with increased radial force for improving attachment and sealing

► **Radiopaque Marks**

Proximal and distal areas for easing the location and position of stent-graft

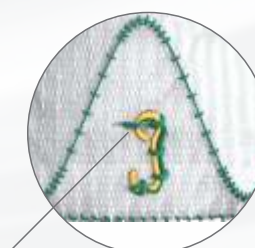

► **Cateter Balão**

Compliant for accommodating stent-graft 12 Fr/45 mm
Supplied separately

Attached by Radial Force
Recommended oversize of 15 to 20% concerning the attaching neck diameter