


## VascoStim Bipolar

Intracardiac electrodes for temporary pacing and electrophysiology

The VascoStim catheters are intravascular catheters for electrophysiological diagnostics and for temporary pacing of the heart, especially in emergencies such as bradyarrhythmia and cardiac conduction system disorder. They can be easily inserted via a peripheral vein.

VascoStim can be positioned quickly, easy and safely and the advantages are good pushability, excellent torque control and positional - stability. VascoStim can be easily user preformed or ordered in traditional shapes.

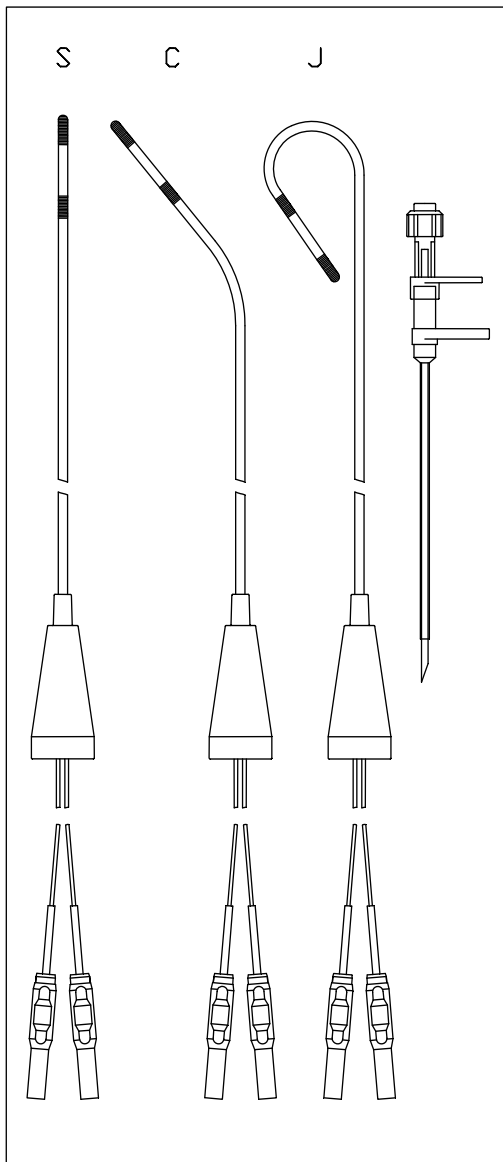

### Technical Data:

Diameter:	5F/1.67mm, 6F/2.0mm
Surface	
Ring/Tip:	15/19mm <sup>2</sup> , 19/25mm <sup>2</sup>
Spacing of Electrodes:	7mm
Material:	Stainless Steel
Eff. length:	114cm
Isolation:	PEBA, radiopaque
Resistance:	< 60 Ohm
<b>Connector:</b>	shrouded 2mm pins

### Attention !

To connect older devices with 2mm connector receptacles please use 2mm adapters (VascoStim TOPAS).

### Package contents:

- 1 Bipolar temporary pacing lead
- 1 Introducer Set

**Packaging Unit:** 5 pcs. / Box

### Ordering Informations:

	Product		
501066	VascoStim	2/5F	C
501065	VascoStim	2/5F	J
501061	VascoStim	2/5F	S
501068	VascoStim	2/6F	C
501067	VascoStim	2/6F	J
501062	VascoStim	2/6F	S

This product is double packaged and sterile.




## Temporary Leads / Electrophysiology

### VascoStim TOPAS 2mm Adapter System

VascoStim TOPAS is an adapter system for temporary pacemakers and EP-Lab systems, so that catheters with 2mm touch protected plugs can be used.

The components of the system allow an easy and safe adaptation of catheters with touch protected 2mm pin connectors to devices with old 2mm connectors (C), as well as a permanent modification of old 2mm pin hole devices to the new standard by means of screwfixable plugs (A).

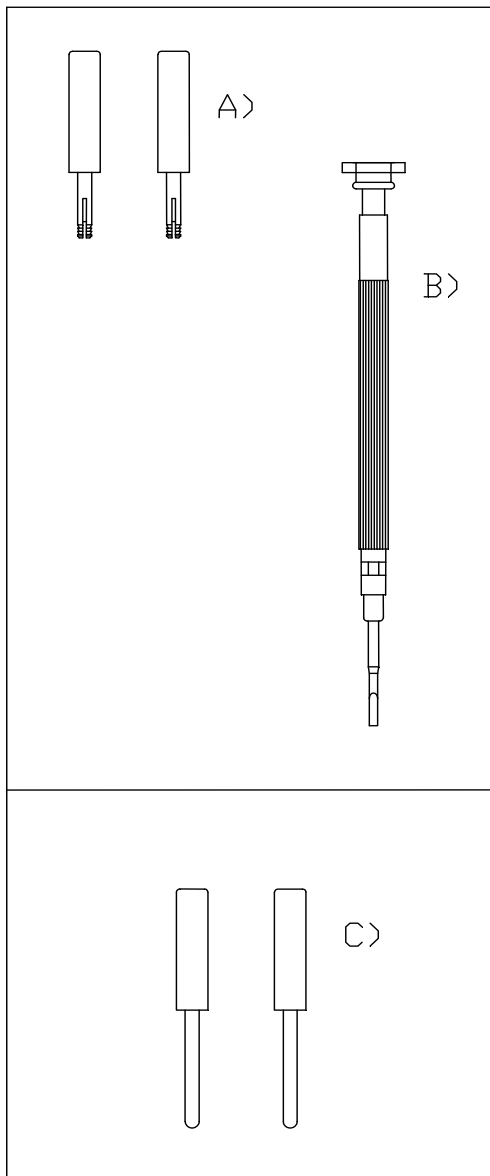

#### Ordering Information:

- 377981 VascoStim TOPAS 5 Set**  
1 Set consisting of:  
- 1 pair Connectors (black/red, C)  
5 Sets per Box
- 509393 VascoStim TOPAS Pacemaker-Set**  
1 Set consisting of:  
- 2 pairs Screwfix plug (black/red, A)  
- 1 pc. Screwdriver (B)  
5 Sets per Box
- 509394 VascoStim TOPAS EP-Set**  
1 Set consisting of:  
- 6 pairs Screwfix plug (black/red, A)  
- 1 pc. Screwdriver (B)  
1 Set per Box
- 509395 VascoStim TOPAS Connector-Set**  
1 Set consisting of:  
- 12 pairs Connectors (black/red, C)  
1 Set per Box
- 509396 VascoStim TOPAS Universal-Set**  
1 Set consisting of:  
- 12 pairs Screwfix plugs (black/red, A)  
- 1 pc. Screwdriver (B)  
- 12 pairs Connectors (black/red, C)  
1 Set per Box