

TEMPORARY BIPOLAR PACING CATHETER

TEMPORARY BIPOLAR PACING CATHETER STANDARD OR SOFT TIP – TORQUE DIRECTED

- intravenous electrode for pre- and postoperative and emergency stimulation
- interim support during pacemaker exchange
- intracardial monitoring
- soft tip version with flexible distal end for atraumatic placement in delicate vessel structures

TEMPORARY BIPOLAR BALLOON PACING CATHETER FLOW DIRECTED


- intravenous placement without fluoroscopic monitoring
- allows for placement at bedside
- for pre- and postoperative and emergency stimulation
- interim support during pacemaker exchange
- intracardial monitoring

FEATURES

- unique design: particularly soft, flexible tip with solid catheter body for excellent torqueability
- antithrombogenic polyurethane catheter body
- available with shrouded (touch protected) pins
- radiopaque
- 10 cm depth markings facilitate precise and safe positioning




shrouded connectors

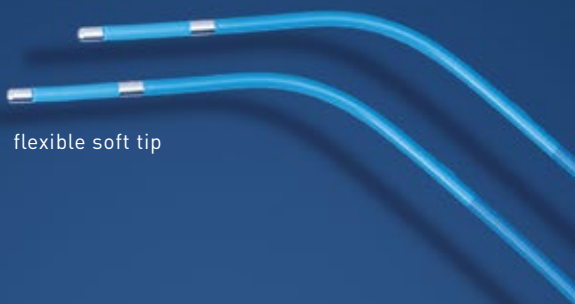

flexible soft tip



depth markings

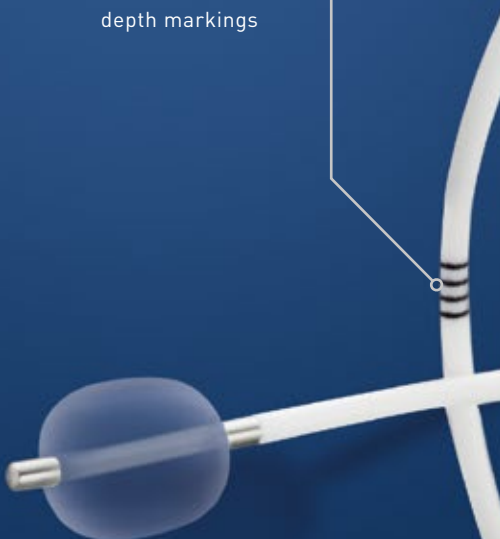

flow directed bipolar
pacing catheter

TEMPORARY BIPOLAR PACING CATHETER

BIPOLAR ELECTRODE FOR TEMPORARY INTRACARDIAC PACING

FIELDS OF APPLICATION:

- Cardiology
- Radiology
- Anaesthesia
- Emergency

PACKAGING:

double packed, sterile

MATERIAL:

Catheter tube: PUR, radiopaque

Electrodes: stainless steel

Balloon: Latex

DIMENSIONS:

Usable length: 110 cm

Depth markings: every 10 cm

Length of flexible distal part: 70 mm

Number of electrodes: 2

Electrode distance: 10 mm

TEMPORARY BIPOLAR PACING CATHETER, TORQUE DIRECTED

Pacing Catheter with touch protected Connectors, FEP Introducer Cannula

Ref.No. Standard Tip	French Size	Electrode Curve	Rcmd. Introducer Size	Max. Balloon Diameter	Packing Unit
24.40110000	4 F / 1.3 mm	Cournand	5 F	-	5 each
24.50110000	5 F / 1.7 mm	Cournand	5 F	-	5 each
24.60110000	6 F / 2.0 mm	Cournand	6 F	-	5 each

Pacing Catheter with touch protected Connectors, without Introducer Cannula

24.70100000	7 F / 2.3 mm	Cournand	7 F	-	5 each
-------------	--------------	----------	-----	---	--------

BIPOLAR PACING CATHETER SETS, TORQUE DIRECTED

Pacing Catheter with touch protected Connectors, Introducer Set with Hemostasis Valve:
Dilator 18 cm / Sheath 11 cm, Seldinger Cannula 18G, Guidewire (double ended) .035" x 45 cm

Ref.No. Standard Tip	French Size	Electrode Curve	Rcmd. Introducer Size	Max. Balloon Diameter	Packing Unit
24.40151100	4 F / 1.3 mm	Cournand	5 F	-	3 each
24.50151100	5 F / 1.7 mm	Cournand	5 F	-	3 each
24.60151100	6 F / 2.0 mm	Cournand	6 F	-	3 each
24.70151100	7 F / 2.3 mm	Cournand	7 F	-	3 each

*Optionally available with sterility sleeve

BIPOLAR BALLOON PACING CATHETER, FLOW DIRECTED

Balloon Pacing Catheter with touch protected Connectors, Syringe

Ref.No.	French Size	Electrode Curve	Rcmd. Introducer Size	Max. Balloon Diameter	Packing Unit
024.BP51	5 F / 1,7 mm	Straight	6 F	9 mm	1 each

ADAPTER PLUG FOR STIMULATION ELECTRODES

FIELDS OF APPLICATION:

- Cardiology
- Radiology
- Anaesthesia
- Emergency

MATERIAL:


Pin: Gold-plated brass
Hub: Polyamide

PACKAGING:

non sterile

ADAPTER PLUG

For safe and easy connection to old type external pacemakers with non-shrouded receptacles.

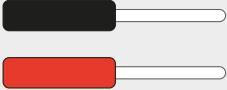Ref.No.		Suitable for Plug Size	Length	Pin Length	Packing Unit
20.24001		ø 2 mm	36,5 mm	18,0 mm	1 pair / bag