




Improved
crossability



Excellent
deliverability

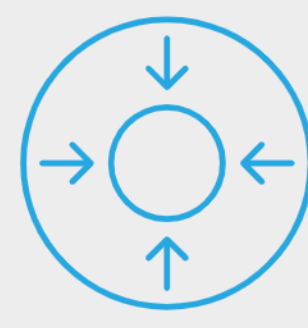

Low profile, wide
range of sizes



Technical data/
ordering info

Teleflex[™]
Empowering the future of healthcare

Passeo[™]-35 Xeo[™]

PTA Balloon Catheter

Low profile PTA balloon with excellent deliverability
in a wide range of sizes.






Passeo™-35 Xeo™ Balloon

Low profile PTA balloon with excellent deliverability in a wide range of sizes.

Improved crossability

Up to 0.16 mm lower tip entry profile compared to leading competitors.¹

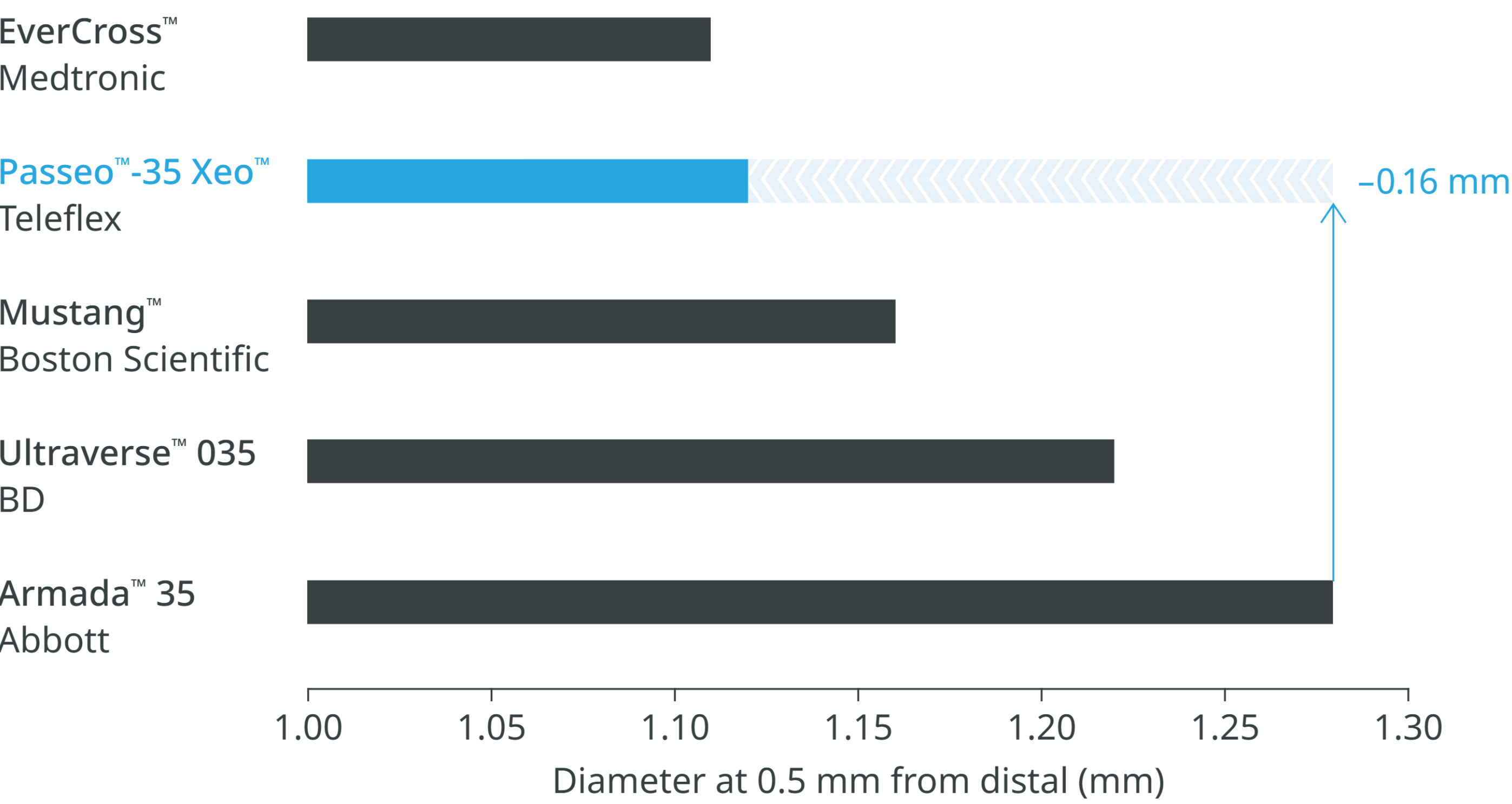


Advanced tip profile

The low tip entry profile and smooth tapered tip of Passeo™-35 Xeo™ Balloon facilitates crossing.



Lowest crossing profile


Up to 0.22 mm lower crossing profile compared to leading competitors.¹

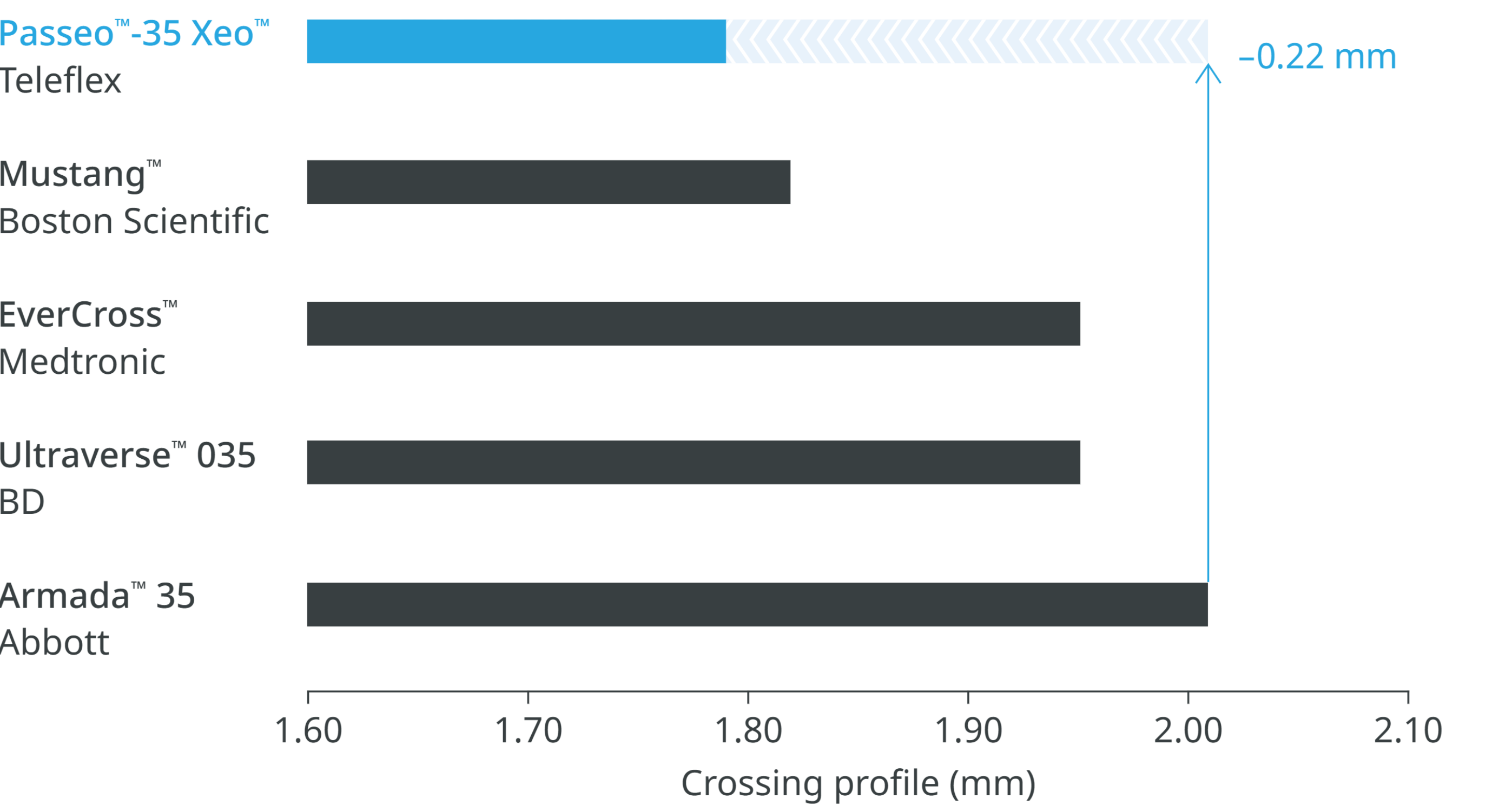



Excellent deliverability

The Passeo™-35 Xeo™ Balloon Catheter offers excellent trackability and pushability compared to competitive 0.035" PTA balloon catheters.²

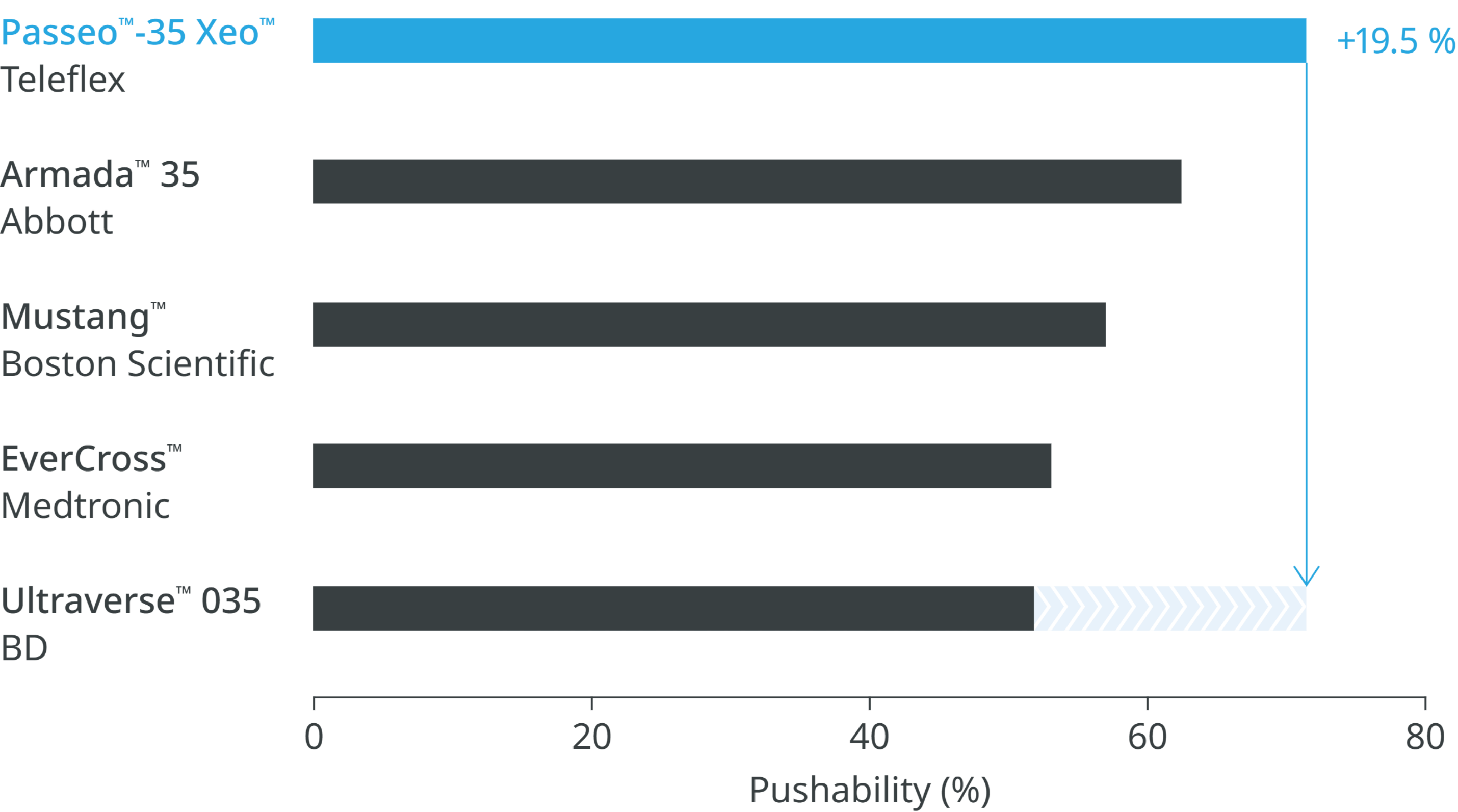


Low profile, wide range of sizes

Passeo™-35 Xeo™ Balloon features a low profile design with a full 5F compatibility up to 7.0 mm balloon diameter and 250 mm balloon length, 6F up to 10.0 mm and 7F up to 12.0 mm diameter.





Passeo™ -35 Xeo™ Balloon
Low profile PTA balloon designed to cross

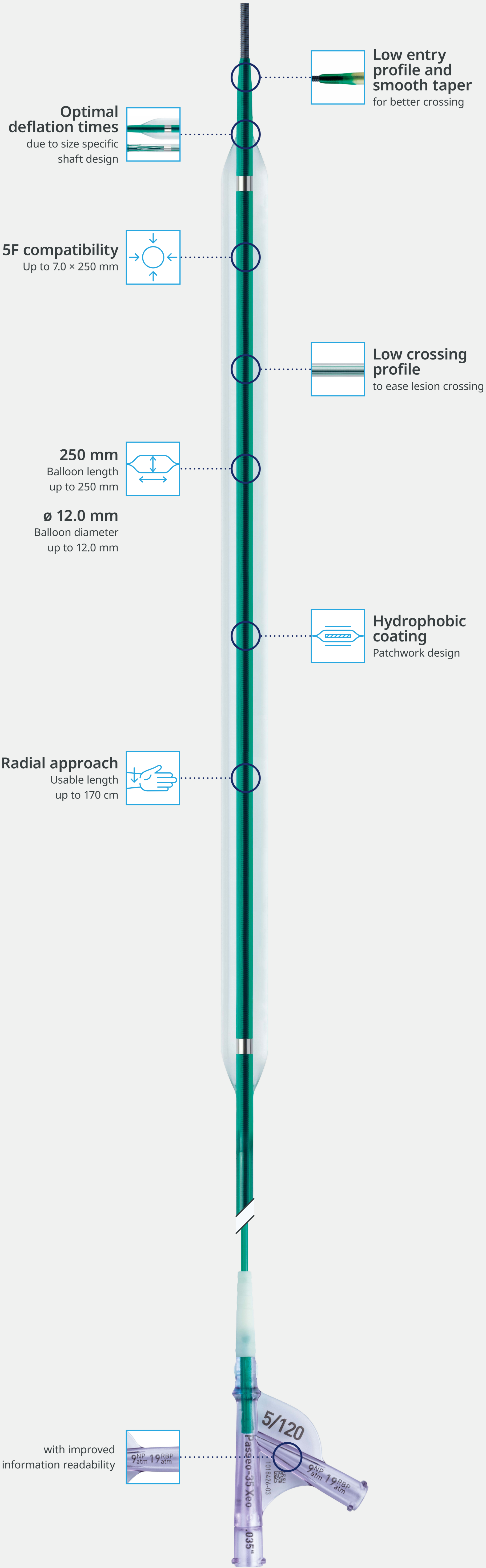

Passeo™-35 Xeo™ Balloon

The Passeo™-35 Xeo™ Balloon is indicated to dilate stenosis in the iliac*, femoral, popliteal and infra-popliteal arteries and for the treatment of obstructive lesions of native or synthetic arteriovenous dialysis fistulae. Passeo™-35 Xeo™ Balloon is also recommended for post-dilation of balloon expandable and self-expanding stents in the peripheral vasculature.**

Vascular Intervention
Peripheral



Technical data	
BALLOON CATHETER	
Catheter type	OTW
Recommended guide wire	0.035"
Tip	Low entry profile, colored
Balloon material	SCP (Semi-Crystalline Polymer), controlled compliance
Balloon coating	Hydrophobic patchwork coating
Balloon markers	2 swaged markers
Sizes	ø 3.0–12.0 mm; L: 20–250 mm
Shaft	5.1-5.4F, dual-lumen, hydrophobic coating
Usable length	90, 130 and 170 cm
Guide wire lumen	Hydrophobic coating

Ordering Information

BALLOON Ø		NP ¹	RBP ²	CATHETER LENGTH 80 CM, BALLOON LENGTH				
				20 mm	40 mm	60 mm	80 mm	100 mm
5F	3.0 mm	9	21	428777	428786	428795	428804	-
	4.0 mm	9	20	428778	428787	428796	428805	428814
	5.0 mm	9	19	428779	428788	428797	428806	428815
	6.0 mm	9	16	428780	428789	428798	428807	428816
	7.0 mm	9	15	428781	428790	428799	428808	428817
6F	8.0 mm	9	14	428782	428791	428800	428809	428818
	9.0 mm	9	12	428783	428792	428801	428810	428819
	10.0 mm	8	11	428784	428793	428802	428811	428820
7F	12.0 mm	8	10	-	428794	428803	428812	428821

BALLOON Ø		NP ¹	RBP ²	CATHETER LENGTH 80 CM, BALLOON LENGTH				
				120 mm	150 mm	170 mm	200 mm	250 mm
5F	3.0 mm	9	21	-	-	-	-	-
	4.0 mm	9	20	428823	-	-	-	-
	5.0 mm	9	19	428824	428833	-	428843	428848
	6.0 mm	9	16	428825	428834	-	428844	428849
	7.0 mm	9	15	-	-	-	-	-
6F	8.0 mm	9	14	-	-	-	-	-
	9.0 mm	9	12	-	-	-	-	-
	10.0 mm	8	11	-	-	-	-	-
7F	12.0 mm	8	10	-	-	-	-	-

BALLOON Ø		NP ¹	RBP ²	CATHETER LENGTH 130 CM, BALLOON LENGTH				
				20 mm	40 mm	60 mm	80 mm	100 mm
5F	3.0 mm	9	21	428851	428860	428869	428878	428887
	4.0 mm	9	20	428852	428861	428870	428879	428888
	5.0 mm	9	19	428853	428862	428871	428880	428889
	6.0 mm	9	16	428854	428863	428872	428881	428890
	7.0 mm	9	15	428855	428864	428873	428882	428891
6F	8.0 mm	9	14	428856	428865	428874	428883	428892
	9.0 mm	9	12	428857	428866	428875	428884	428893
	10.0 mm	8	11	428858	428867	428876	428885	428894
7F	12.0 mm	8	10	-	428868	428877	428886	428895

BALLOON Ø		NP ¹	RBP ²	CATHETER LENGTH 130 CM, BALLOON LENGTH				
				120 mm	150 mm	170 mm	200 mm	250 mm
5F	3.0 mm	9	21	428896	428905	428910	428915	428920
	4.0 mm	9	20	428897	428906	428911	428916	428921
	5.0 mm	9	19	428898	428907	428912	428917	428922
	6.0 mm	9	16	428899	428908	428913	428918	428923
	7.0 mm	9	15	428900	428909	428914	428919	428924
6F	8.0 mm	9	14	428901	-	-	-	-
	9.0 mm	9	12	428902	-	-	-	-
	10.0 mm	8	11	428903	-	-	-	-
7F	12.0 mm	8	10	428904	-	-	-	-

BALLOON Ø		NP ¹	RBP ²	CATHETER LENGTH 130 CM, BALLOON LENGTH				
				20 mm	40 mm	60 mm	80 mm	100 mm
5F	4.0 mm	9	20	-	428935	428944	428953	-
	5.0 mm	9	19	428927	428936	428945	428954	-
	6.0 mm	9	16	428928	428937	428946	428955	-
	7.0 mm	9	15	428929	428938	428947	428956	-
	8.0 mm	9	14	428930	428939	428948	428957	-
6F	9.0 mm	9	12	-	428940	428949	428958	-
	10.0 mm	8	11	-	428941	428950	428959	-

BALLOON Ø		NP ¹	RBP ²	CATHETER LENGTH 130 CM, BALLOON LENGTH				
				120 mm	150 mm	170 mm	200 mm	250 mm
5F	4.0 mm	9	20	428971	428980	-	428990	428995
	5.0 mm	9	19	428972	428981	-	428991	428996
	6.0 mm	9	16	428973	428982	-	428992	428997
	7.0 mm	9	15	-	-	-	-	-
6F	8.0 mm	9	14	-	-	-	-	-
	9.0 mm	9	12	-	-	-	-	-
	10.0 mm	8	11	-	-	-	-	-

¹ NP=Nominal Pressure and RBP=Rated Burst Pressure in atm (1 atm=1.013 bar)
^{*} Note for Australia: Passeo™-35 Xeo™ Balloon not approved by Therapeutic Goods Administration for use in the common iliac arteries.
^{**} Indication as per IFU.