


Exeter Retrieval Snare

fingerspitzengefühl in fishing for foreign objects

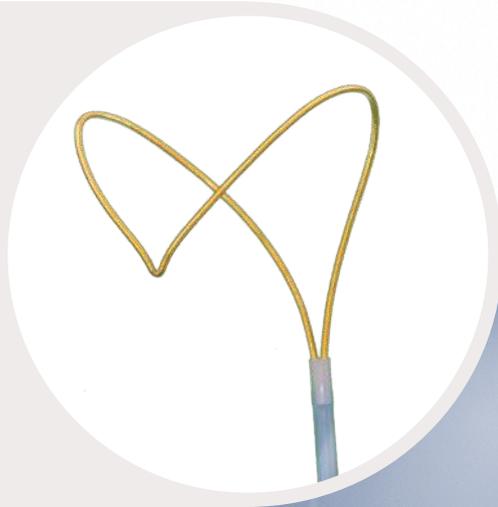

Exeter Retrieval Snare

- 1 Atraumatic Nitinol loop with gold-plated tungsten for highest visibility.
- 2 Pre-shaped MP Tip eases navigation.
- 3 Shapable introducer to ease navigation and gain performances.
- 4 AndraTec Supportive sleeve avoids kinking when manipulating.
- 5 Luer Connector.
- 6 Extra XXL torque device to ease manipulation
- 7 One piece kink-resistant wire enhances 1:1 torque for better manipulation.
- 8 Radiopaque platinum inner Marker Ring ensures nearest positioning of introducer.



The **Exeter Snare** and **Exeter Snare Micro** have been especially designed to meet the high requirements during extraction.

The **Exeter Retriever Technology** is your first choice safety device when it comes to successful retrieval.

Order codes

Product Code	⌀ (mm)		↓ (cm)	⊔ (cm)	⊔ (F)
ES-05	5mm	(3-5mm vessel)	125	110	4
ES-10	10mm	(6-10mm vessel)	125	110	4
ES-15	15mm	(11-15mm vessel)	125	110	4
ES-20	20mm	(16-20mm vessel)	125	110	4
ES-25	25mm	(21-25mm vessel)	125	110	5
ES-30	30mm	(26-30mm vessel)	125	110	5
ES-35	35mm	(31-35mm vessel)	125	110	5

Product Code	⌀ (mm)		↓ (cm)	⊔ (cm)	⊔ (F)
ES-25/65	25mm	(21-25mm vessel)	65	50	5
ES-25/150	25mm	(21-25mm vessel)	150	135	5

Product Code	⌀ (mm)		↓ (cm)	⊔ (cm)	⊔ (F)
ESM-2	2mm	(1-2mm vessel)	175	150	2-3-2
ESM-4	4mm	(3-4mm vessel)	175	150	2-3-2
ESM-7	7mm	(5-7mm vessel)	175	150	2-3-2



Andratec GmbH
 Simmener Straße. 70, D-56075 Koblenz, Germany
Office & Warehouse
 Im Deutschen Kaiser,
 Kastorstraße 3, 56068 Koblenz, Germany
 ☎ +49(0) 261 100 49 653
 📠 +49(0) 261 100 49 654
 📧 info@andratec.com
 🌐 www.andratec.com

